

CD3E & CD3D Heterodimer protein-coupled magnetic beads

Cat. No. CD3-003M **Lot. No.** (See product label)

SPECIFICATION

Product Overview

The biotinylated CD3E & CD3D protein was conjugated to streptavidin magnetic beads. This pre-coupled magnetic bead product can capture the anti-CD3E & CD3D antibody from various assay systems. The beads are in uniform size, narrow size distribution with large surface area and unique surface coating, which can help you get the best performance and highly reproducible results. This Human CD3E & CD3D coupled magnetic beads will bring great convenience with minimum non-specific binding and developed protocols. This ready-to-use product could greatly save your time and hassle.

Species Human

Beads Size 2mg

Particle size 2 μ m

Beads Surface hydrophilic

Coupled amount of protein 30 μ g CD3E & CD3D /mg Beads


Capacity >200 nmol antibody/ mg beads

Formulation Lyophilized from 0.22 μ m filtered solution in PBS, 0.05% Tween-20, pH7.4, with 10% trehalose


 Tel: 1-631-559-9269 1-516-512-3133

 Email: info@creative-biomart.com  Fax: 1-631-938-8127

 45-1 Ramsey Road, Shirley, NY 11967, USA

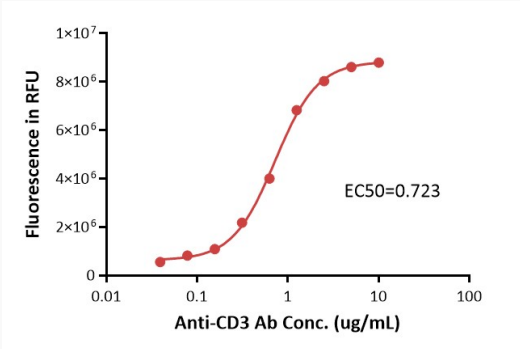Reconstitution	2 mL ultrapure water (1mg beads/mL)
Application	This product is intended for immunocapture, biopanning and flow cytometry.
Storage	Upon receipt, please store the Beads at -20°C for 1 year in lyophilized state. Once the Beads reconstitution, please use it immediately. Do not to freeze thaw the Beads after reconstitution.
Assay Principles	The conjugation was achieved by means of the binding between streptavidin and biotin. Streptavidin (SA) has an extraordinarily high affinity for biotin with a dissociation constant (Kd) on the order of 10-14 mol/L. Thus, the binding of streptavidin and biotin is irreversible. Our Human CD3E & CD3D-coupled magnetic beads could capture anything binding to CD3E & CD3D, and make the following testing easy, such as immunocapture, biopanning and flow cytometry.
Application Method	a) Reconstitute the Beads following the COA. Wash and re-suspended the beads to a certain concentration by adding your dilution buffer. b) Add the prepared beads to your samples. c) Beads can be separated from your samples afterwards using a magnetic plate.
	 Add CD3 pre-coupling magnetic beads into each well at a concentration of 1 mg/ml.

0.1 mg beads was added to each well. Then add increasing concentrations of anti-CD3 antibody to different wells. PE-labeled anti-human IgG Fc antibody was used for detection.

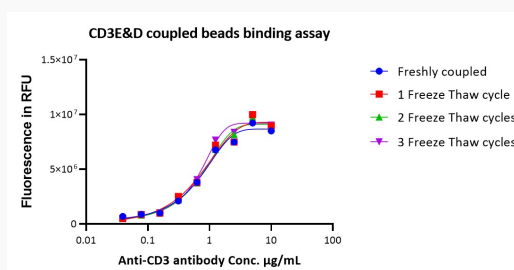

Add CD3 pre-coupling magnetic beads into each well at a concentration of 1 mg/ml. 0.1 mg beads was added to each well. Then add increasing concentrations of anti-CD3 antibody to different wells. PE-labeled anti-human IgG Fc antibody was used for detection.