

Whole Blood CD15 protein-coupled magnetic MicroBeads

Cat. No. CD15-7453M **Lot. No.** (See product label)

SPECIFICATION

Product Overview

Species Human

Capacity for 40 mL whole blood


Background

StraightFrom Whole Blood CD15 MicroBeads recognize the carbohydrate structure 3-fucosyl-N-acetylglucosamine (3-FAL), also designated Lewis X. The CD15 antigen is expressed on neutrophils and eosinophils but not on basophils and lymphocytes.

Application

StraightFrom Whole Blood MicroBeads allow the fast isolation of CD15+ cells directly from whole blood or bone marrow, thus minimizing hands-on time and maximizing yield of target cells. StraightFrom Whole Blood MicroBeads in combination with the autoMACS Separators allow standardization of the cell separation procedure and safe handling of hazardous samples. The purified CD15+ cells are well-suited for further flow cytometric, functional, or molecular analysis including lineage-specific chimerism analysis after allogeneic stem cell transplantation.

Analysis

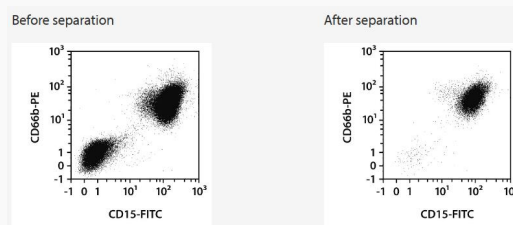

Separation of a whole blood sample using the StraightFrom Whole Blood CD15

 Tel: 1-631-559-9269 1-516-512-3133

 Email: info@creative-biomart.com  Fax: 1-631-938-8127

 45-1 Ramsey Road, Shirley, NY 11967, USA

MicroBeads and the MultiMACS Cell24 Separator Plus with the Single-Column Adapter and Whole Blood Columns without washing step after labeling. Cells were fluorescently stained with CD15-FITC, CD66b-PE, as well as CD45-VioBlue and analyzed by flow cytometry using the MACSQuant Analyzer. Cells were triggered via CD45-VioBlue, cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.

Tel: 1-631-559-9269 1-516-512-3133

Email: info@creative-biomart.com Fax: 1-631-938-8127

45-1 Ramsey Road, Shirley, NY 11967, USA