

Whole Blood REAlease CD8 protein-coupled magnetic MicroBeads

Cat. No. CD8-6428M **Lot. No.** (See product label)

SPECIFICATION

Product Overview

Species Human

Capacity for 40 mL whole blood


Background

StraightFrom Whole Blood REAlease CD8 MicroBead Kit was developed for rapid positive selection of CD8+ cells directly from whole blood. No sample preparation is required, including density gradient centrifugation or erythrocyte lysis. The purified target cells are bead-free, practical for consecutive magnetic isolations, and can be label-free so that epitopes are available for any downstream application.

Application

StraightFrom Whole Blood REAlease CD8 MicroBead Kit allows the fast isolation of CD8+ cells directly from blood thus minimizing hands-on time and maximizing yields of target cells. The purified CD8+ cells are well-suited for further flow cytometric, functional, or molecular analysis. The isolated CD8+ cell subset is label-free.

Analysis

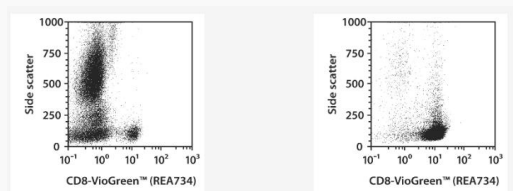

CD8+ cells were isolated from whole blood using the StraightFrom Whole Blood

 Tel: 1-631-559-9269 1-516-512-3133

 Email: info@creative-biomart.com  Fax: 1-631-938-8127

 45-1 Ramsey Road, Shirley, NY 11967, USA

REAl ease CD8 MicroBead Kit, human and Whole Blood Columns. Cells were fluorescently stained with the REAfinity Antibody CD8-FITC and analyzed by flow cytometry using the MACSQuant Analyzer X. Cell debris and dead cells were excluded from the analysis based on scatter signals and propidium iodide fluorescence.

 Tel: 1-631-559-9269 1-516-512-3133

 Email: info@creative-biomart.com  Fax: 1-631-938-8127

 45-1 Ramsey Road, Shirley, NY 11967, USA