

Ubiquitin Magnetic CoIP Agarose

Cat. No. Ubiquitin-04A **Lot. No.** (See product label)

SPECIFICATION

Product Overview	<p>Immunoprecipitation of ubiquitin and ubiquitinated proteins with anti-Ubiquitin Nanobody conjugated to magnetic agarose beads</p> <ul style="list-style-type: none"> • Fast, reliable & efficient one-step immunoprecipitation • Ready-to-use • No heavy & light antibody chains • Stable under harsh washing conditions • Suitable for downstream mass spec analysis • Works in samples from: mammals, plants, yeast, etc.
Specificity/Target	<p>Binds to monomeric ubiquitin, ubiquitin chains, and ubiquitinated proteins. Ubiquitin chains may be linked via lysines; linkages via lysines such as K11, K48 and K63 are compatible with binding by the Ubiquitin-Trap (further linkage types also likely to be compatible, but not tested).</p>
Tested Reactivity	Human, Mouse, Hamster, Dog, Spinach, Yeast
Conjugate	Magnetic agarose beads; ~40 um (cross-linked 6% magnetic agarose beads)
Elution buffer	2xSDS-sample buffer (Lamml)
Wash buffer compatibility	2M NaCl, 5 mM DTT, 0 mM β -mercaptoethanol, 5 mM TCEP, 2% NP40, 2% Triton X-100, 0% SDS, 2-3 M Urea
Type	Nanobody

 Tel: 1-631-559-9269 1-516-512-3133

 Email: info@creative-biomart.com  Fax: 1-631-938-8127

 45-1 Ramsey Road, Shirley, NY 11967, USA

Affinity (KD)

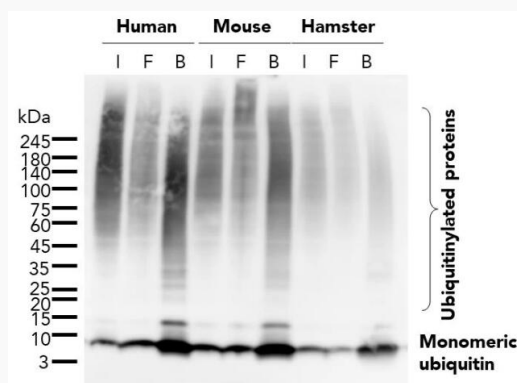
The KD for monomeric ubiquitin is 90 nM. The apparent affinity is even higher for ubiquitin chains due to avidity effects.

Storage Buffer

20% ethanol

Applications

IP, CoIP

Monomeric ubiquitin


Tel: 1-631-559-9269 1-516-512-3133

Email: info@creative-biomart.com Fax: 1-631-938-8127

45-1 Ramsey Road, Shirley, NY 11967, USA